


Leather Color Selection




Option

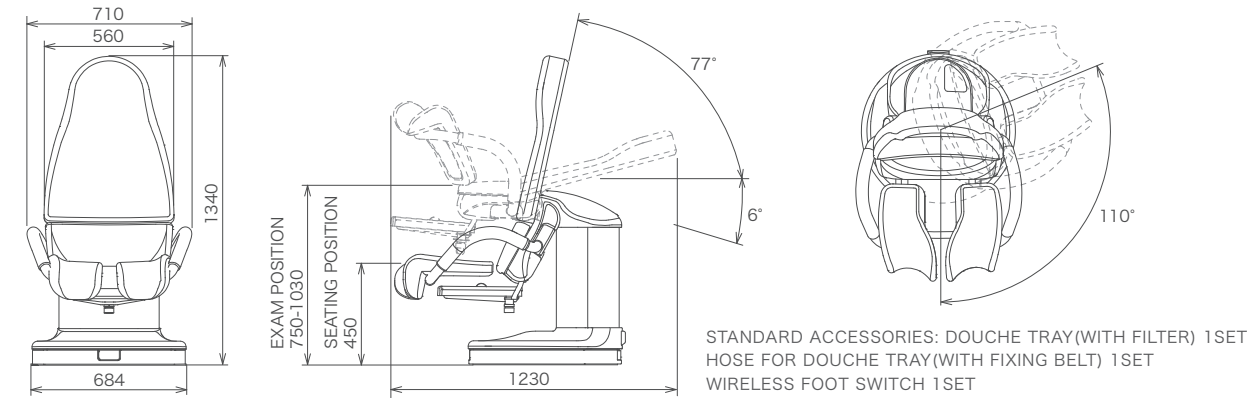

GYNECOLOGICAL UNIT  
DG-100

The DG-100 gyneco logical unit is designed as a low-cost model.  
The thermo-irrigator is optional and is hung on the bar.  
The tabletop is made of stainless steel and can be used as an instrument tray.  
A large drawer and two shelves are available for storage.  
The DG-501 and DG-401 light is optional.

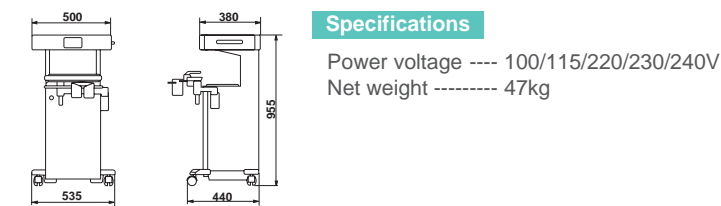

DG-100 with DG-501

GYNECOLOGICAL UNIT  
DG-315

Ease of operation, yet functional  
with high efficiency.




- The reliable electro-hydraulic system eliminates vibration and shock while being operated.
- The table is equipped with a preset and auto return device for speedy examinations and treatments.
- The douche tray (irrigation tray) slides smoothly below the seat.
- The table height can be adjusted between 49.5cm and 99.5cm to provide doctor examining convenience and the patient comfort and confidence.

Color selection  
of the leather






Optional accessories

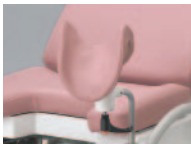

Dimensions

-mm-

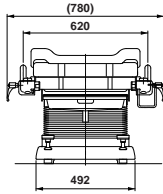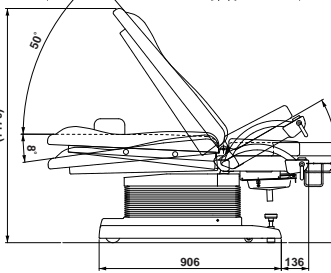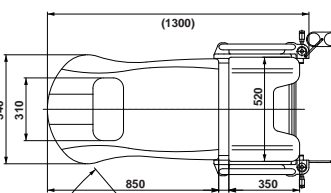

Specifications

- |                                |       |  |
|--------------------------------|-------|--|
| Table lifting mechanism        | ----- | Electro-hydraulic system by Foot control |
| Back support lifting mechanism | --    | Electro-hydraulic system by Foot control |
| Seat tilting mechanism         | ----- | Electro-hydraulic system by Foot control |
| Table lifting                  | ----- | With Preset Device                       |
| Table lowering                 | ----- | With Auto Return                         |
| Power voltage                  | ----- | 115/120/220/230/240 AC                   |
| Control voltage                | ----- | DC-12V                                   |
| Net weight                     | ----- | 160kg                                    |

Standard accessories:

- pillow-1pc, douche tray with wire gauze-1pc, exhaust hose with band-1pc, Foot control-1pc, instrument holder-1pc, cotton container-1pc, forceps jar-1pc, plastic seat cover-1pc

Safety  
Precautions

Confirm the installation conditions, use conditions, and purpose prior to purchase. Be sure to carefully read and follow the attached documents and instructions before use. The dimensions, weight, and working range shown in this brochure include some margins of error. Product specifications and design are subject to change without notice. Colors may look slightly different from those of the actual product due to printed matter.

The manufacturer reserves the right to alter specifications without prior notice.



TAKARA BELMONT CORPORATION  
2-1-1, Higashishinsaibashi,  
Chuo-ku, Osaka, 542-0083, Japan  
TEL : 81-6-6213-5945  
FAX : 81-6-6212-3680

Book no. M-004  
Printed in Japan 2012. 09

OBSTETRIC &  
GYNECOLOGICAL  
EXAMINATION TABLE






# Gynecological Table DG-770

**SIMPLE**  
Roundish shape creates gentleness.

**SPEEDY**  
Quick position change for examination.

**SPACE**  
New rotating mechanism realizes space-saving.

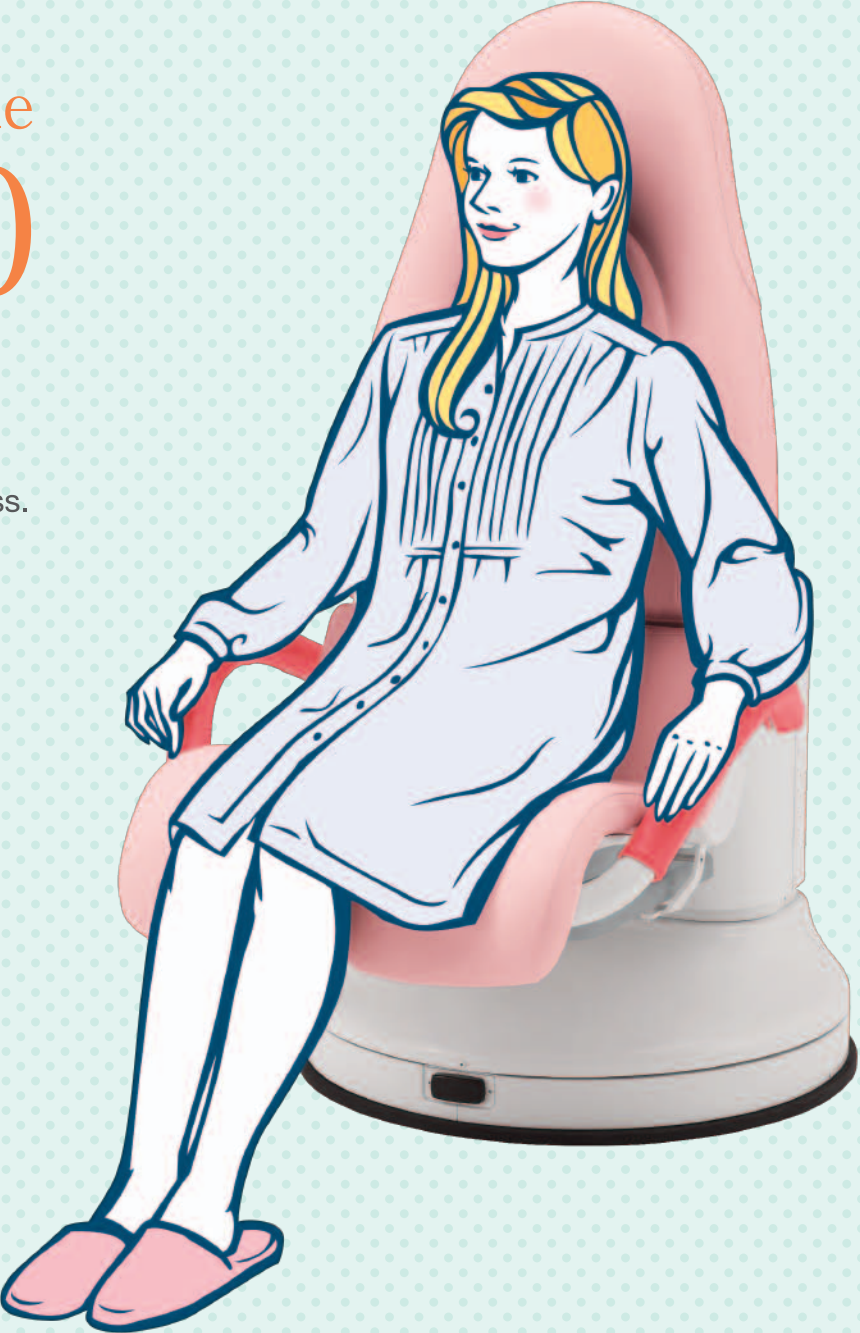

## Pursuing Quality of Time

- To allow for effective use of consultation hours
  - To minimize anxiety of patients
  - To provide time with relaxation
- Gynecological examination chair DG-770 gives consideration to “time and space” on every occasion.



## SIMPLE


**Warmth by shape and softness through the skin**

Warm feeling by looking. Enfolding design gently supports patient's body. Comfortable fit and cushioning properties of thigh rests with a sense of unity allow for gentle support of the hucklebone like enveloping in a warm hug, providing a stable sitting posture.

## SPEEDY

**Smooth movement for quick position change**


Time for change from initial position to examination position is reduced to enable consultation hours to be effectively used. (Compared with our existing products)

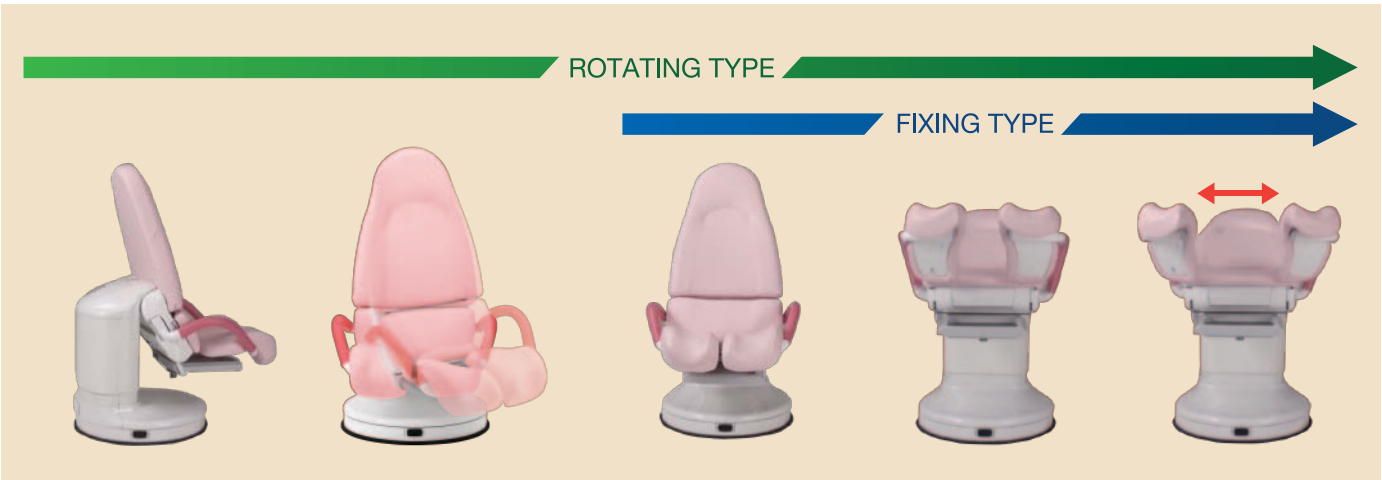

## SPACE


**Compact design with top priority on effective use of consultation rooms**

Reduced rotation radius saves installation space. Minimum installation dimensions are reduced to facilitate installation of rotary type chairs. (Compared with our existing products)

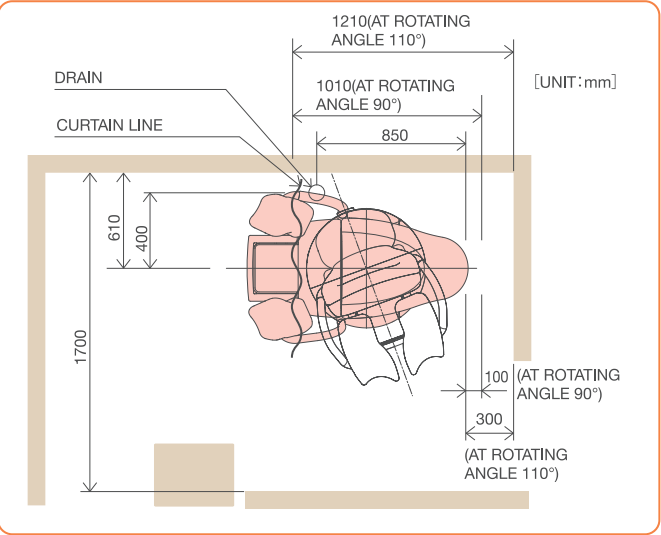

## Useful Features

- Examination area**  
Widened leg opening secures a sufficient examination area. Thigh rests are in shape with cushioning properties maintained and thickness reduced, allowing for smooth transvaginal scan.
- Leg-closed examination position**  
Leg-closed examination position switch allows for position change for examination with legs closed. Individualized leg opening can be operated according to the condition of patients including those after saturation and those with disease in the hip joint.
- Selective leg opening**  
Foot switch allows you to easily open right and left legs individually. You can gradually open one after the other to respond to each patient.
- Wireless foot switch**  
Infrared-based wireless foot switch is employed to require no cords, offering a neat examination room.

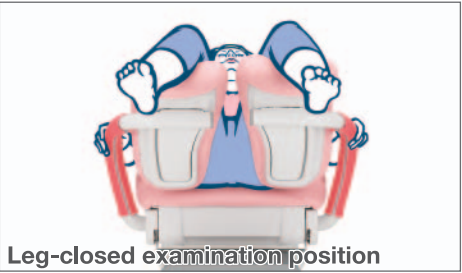

**Plastic douche tray**

Corrosion resistance and water repellent finishing have derived hard-to-stain and easy-to-clean tray. Back-and-forth sliding mechanism allows you to adjust the tray position as necessary.



**Shockless rotation**

Hydraulic pressure control realizes soft starting and soft stopping, improving smooth movement of the examination chair.

**Main unit auxiliary switches**

Main unit is equipped with “P: Preset” and “R: Return” functional switches, allowing for auxiliary operation from the main unit.

